

ANTERION®

Multimodal Imaging Platform
Optimized for the Anterior Segment

HEIDELBERG
ENGINEERING



The multimodal functionality you need

The image quality you expect

Imaging App

- Anterior chamber and angle imaging
- Corneal and scleral imaging
- Visualization of the lens and both surfaces
- Customizable scan patterns
- Peripheral imaging

The expandable design you want (optional)

Cornea App

- Corneal topography
- Corneal tomography
- Pachymetry
- Total corneal power
- Corneal wavefront analysis
- Corneal differential maps
- Progression analysis

Cataract App

- Axial length
- Lens thickness
- Aqueous depth
- Central corneal thickness
- Anterior axial curvature
- Total corneal power
- Total corneal wavefront
- Spheric and toric IOL calculator
- Formulas:
 - SRK/T
 - Haigis
 - Holladay I
 - Hoffer® Q
 - Barrett formulas (Universal II, Toric, True-K, True-K Toric)
 - ULIB import
 - IOL CON import
 - Okulix interface


Metrics App


- Anterior chamber angle assessment
- 360° graphs of angle parameters
- Anterior chamber volume
- Lens vault
- Lens thickness
- Free-hand measurements

Experience confident diagnostics and workflow optimization

- Visually confirm all measurements with exceptionally clear swept-source OCT images.
- Increase efficiency and save space with a customizable, compact platform that meets your individual needs.
- Simplify the patient journey and transform your clinical workflow.
- Discover the speed and security of HEIDELBERG EYE EXPLORER HEYEX 2 image management.


Anterior segment anatomy

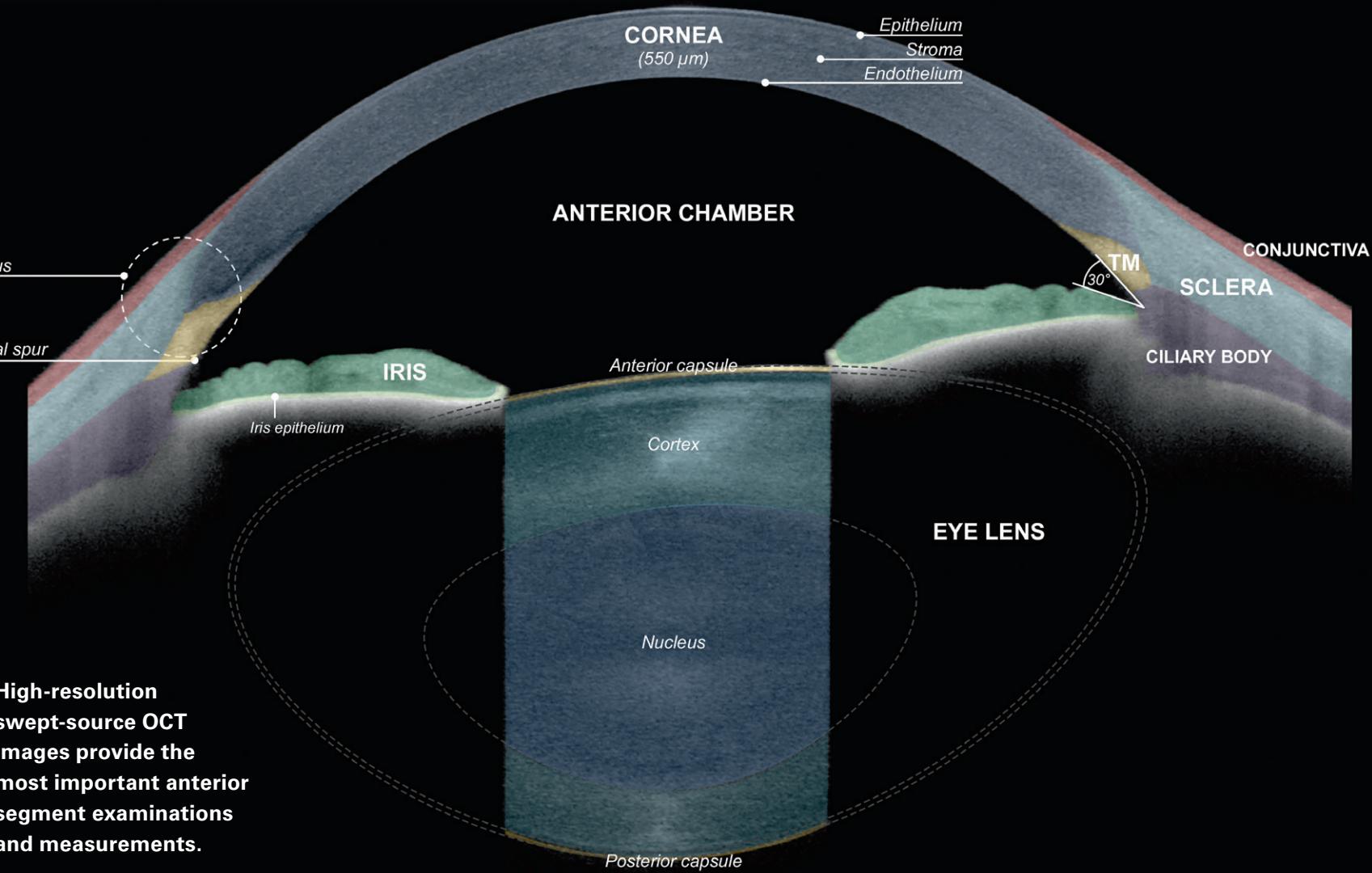

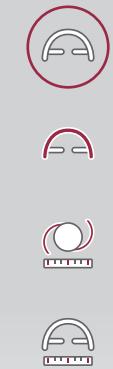

Imaging App - Enhance your clinical confidence


The Imaging App comes as standard. Visualize various anterior segment pathologies or surgical results, such as keratoplasty, implanted IOLs, and corneal rings.

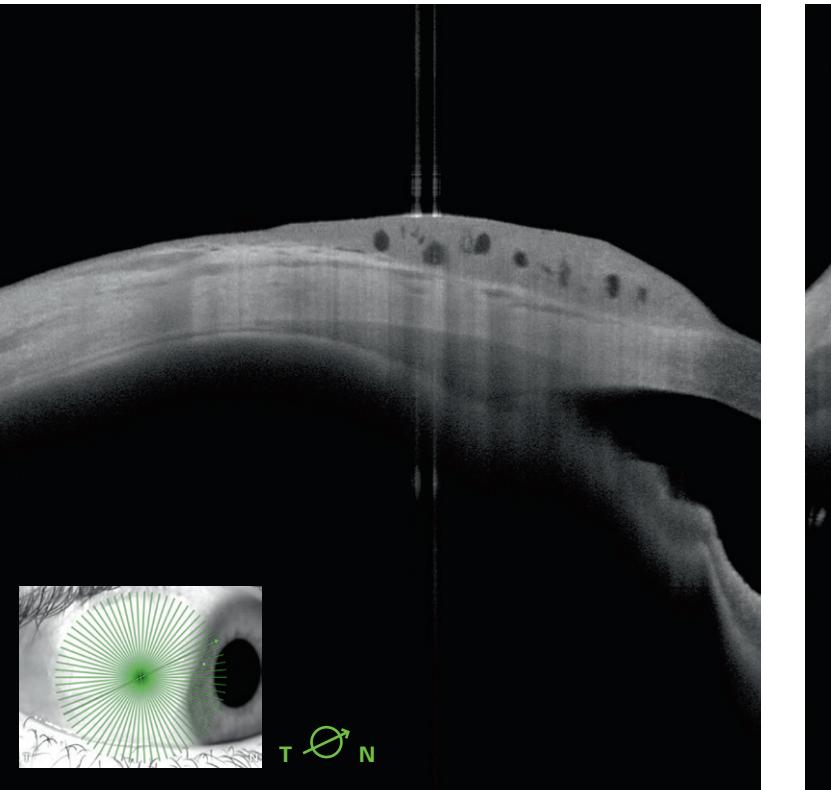


Conjunctival nevus (lateral fixation)

Image courtesy: Sacha Nahon-Esteve, MD, Nice, France

A detailed visualization of the sclera, ciliary body, and rectus muscle assists you in the diagnosis and management of pathologies that impact these anatomical structures.

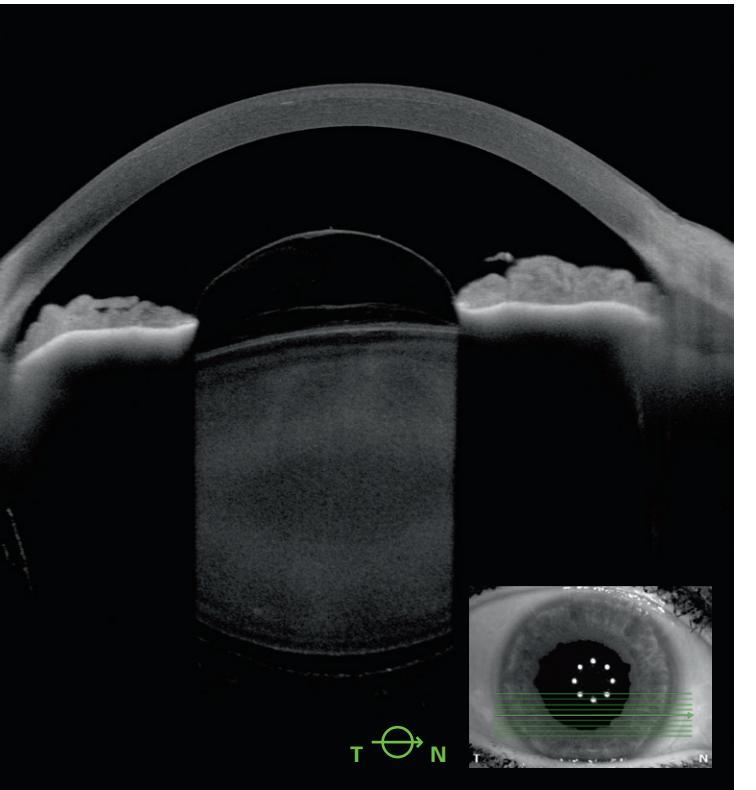



Post trauma, vitreous in the anterior chamber

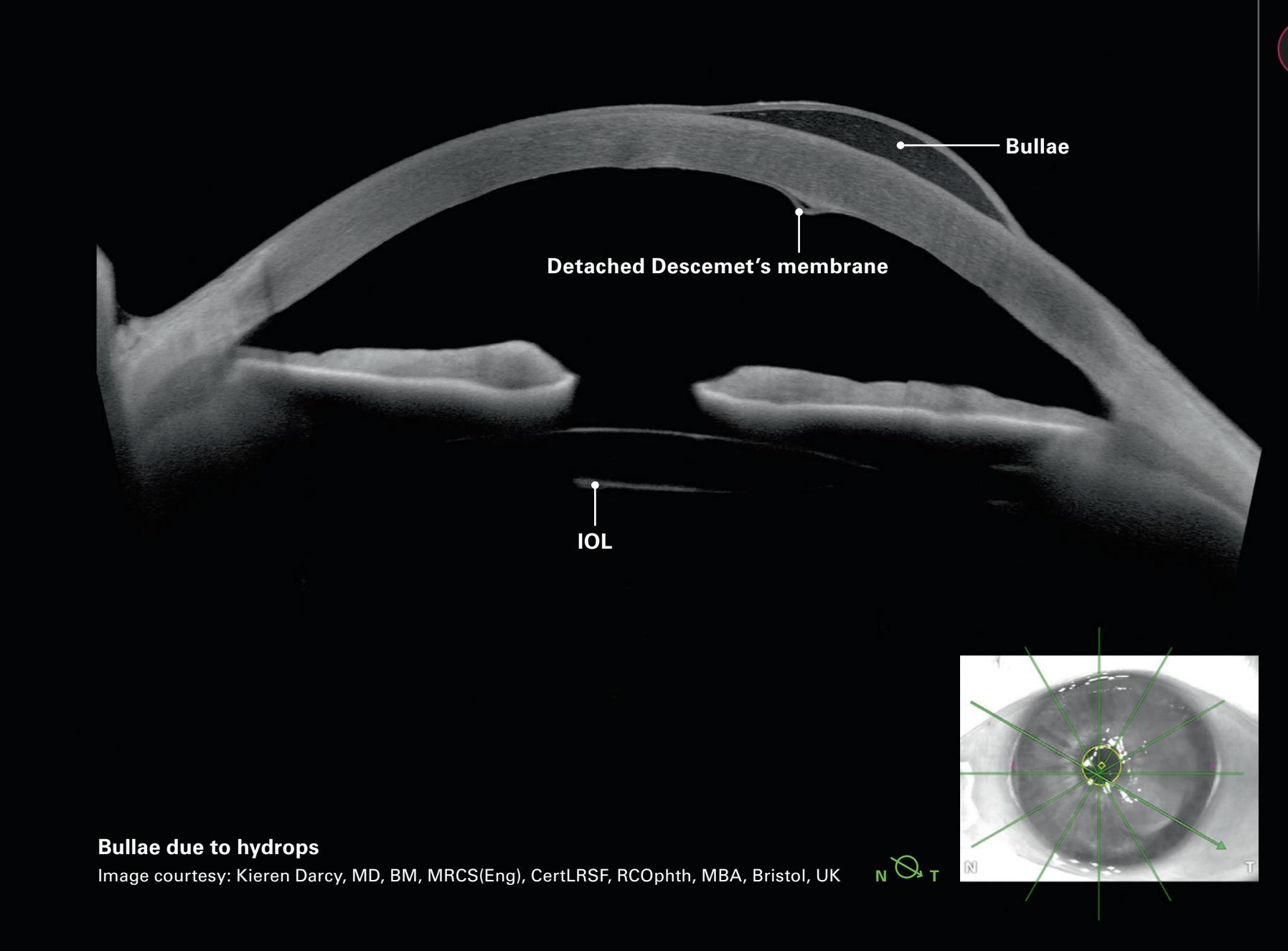



Bullae due to hydrops

Image courtesy: Kieren Darcy, MD, BM, MRCS(Eng), CertLRSF, RCOPhtth, MBA, Bristol, UK



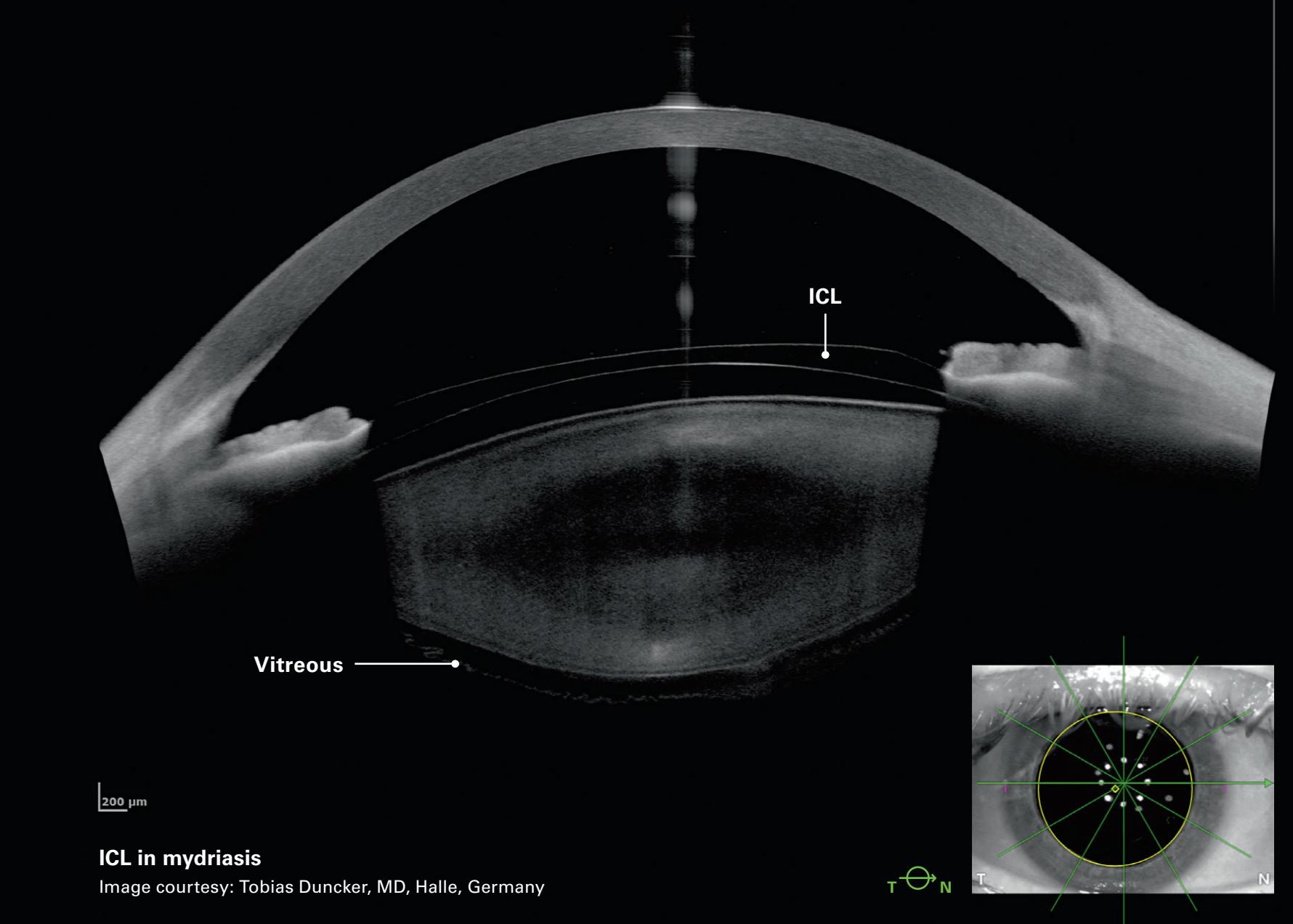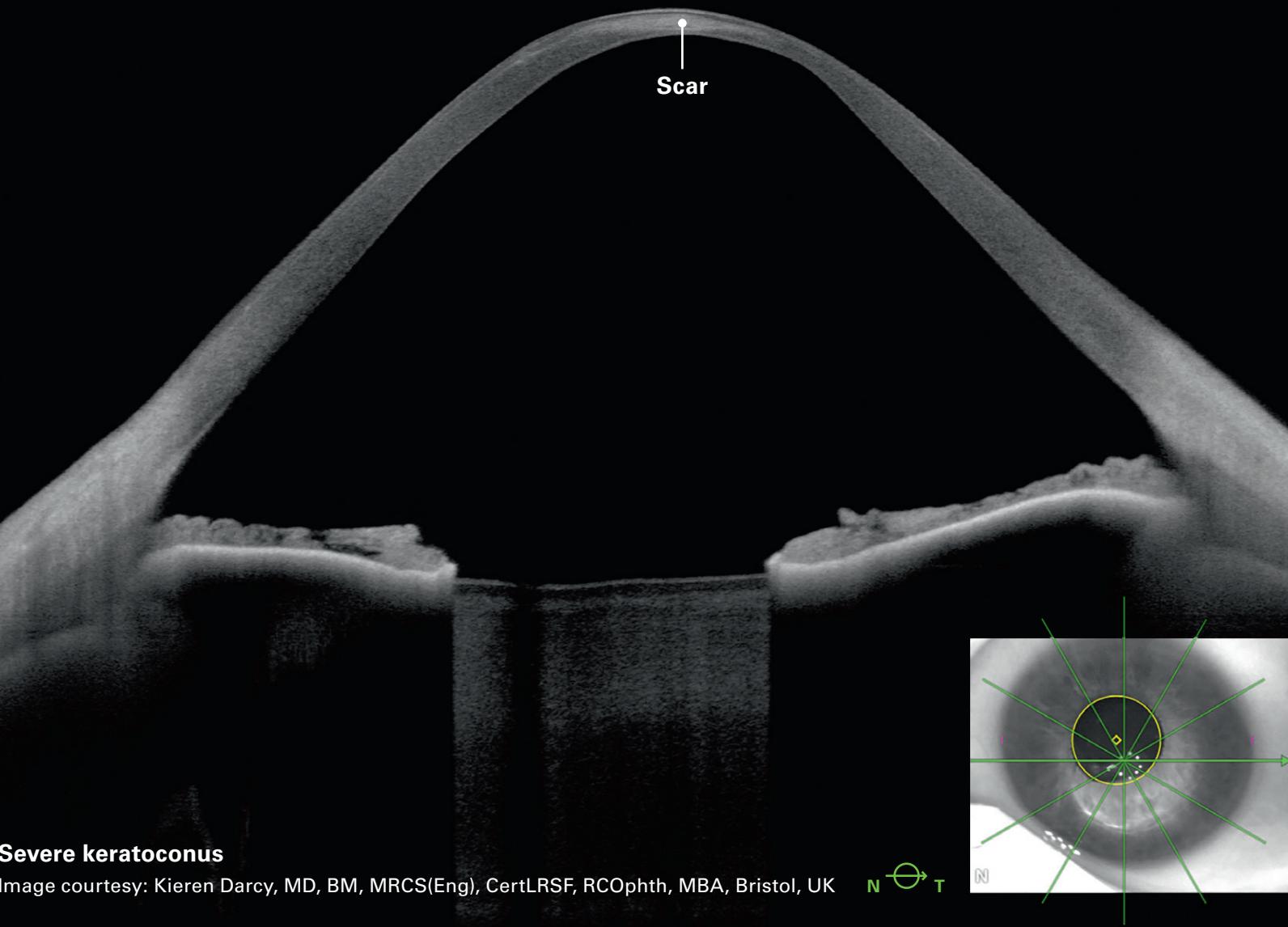

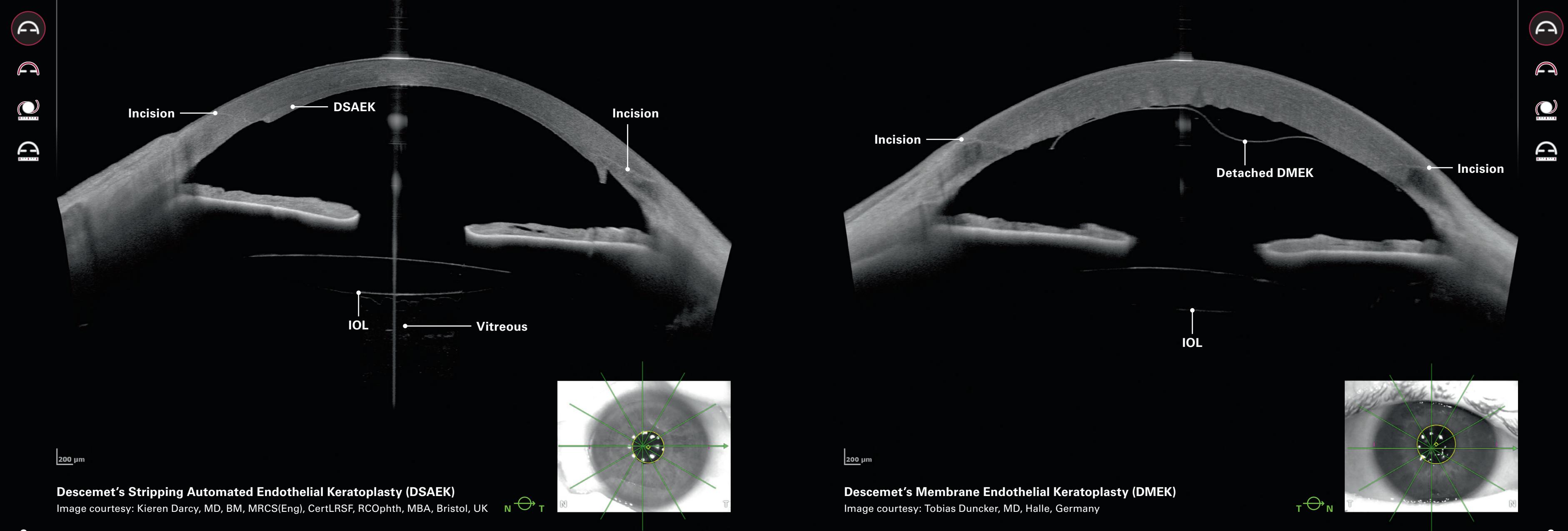


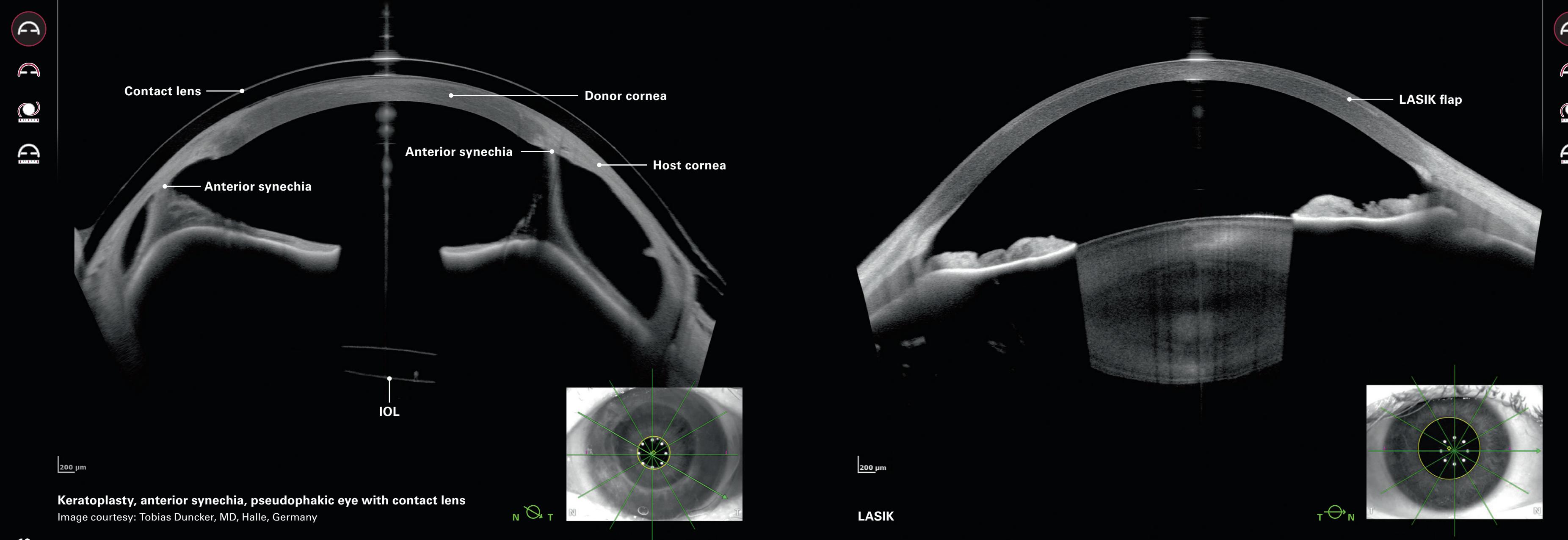




Cornea App – Customized reports enhance your workflow


Anterior axial curvature

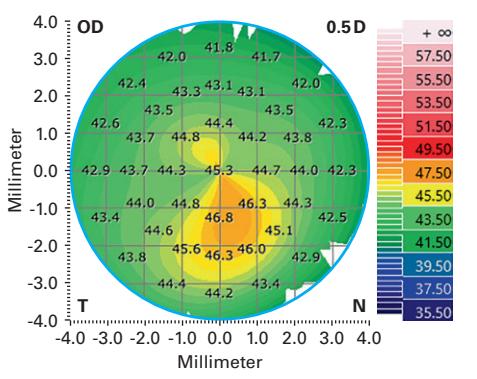


Anterior tangential curvature

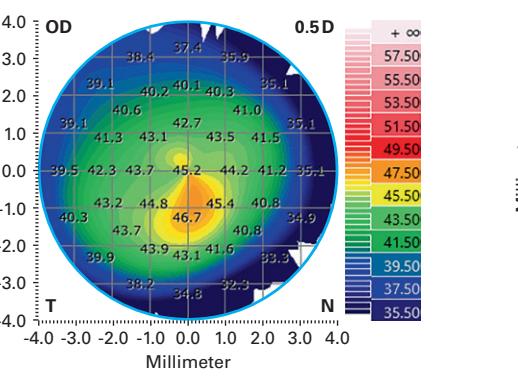


Pachymetry

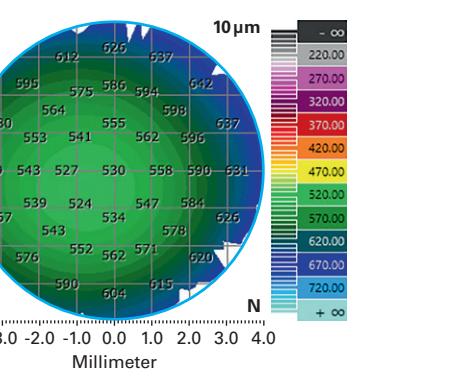


Posterior axial curvature

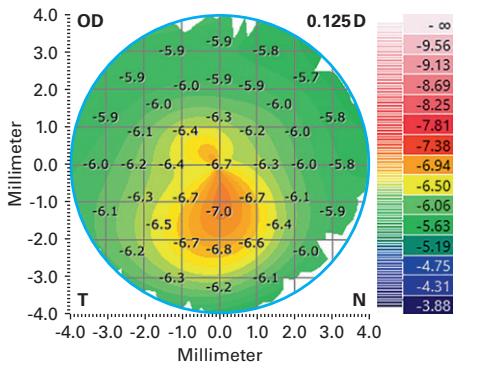

Posterior tangential curvature

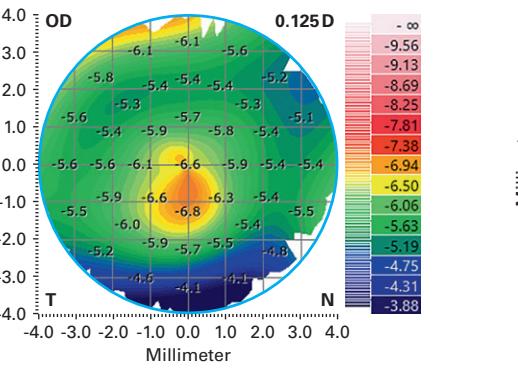

Total corneal power

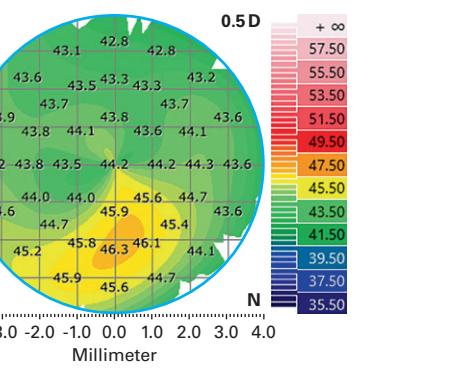


Keratoconus Multiview


Image courtesy: Oliver Findl, MD, MBA, FEBO, Vienna, Austria

Corneal maps correlate with camera and OCT images


Total corneal power map

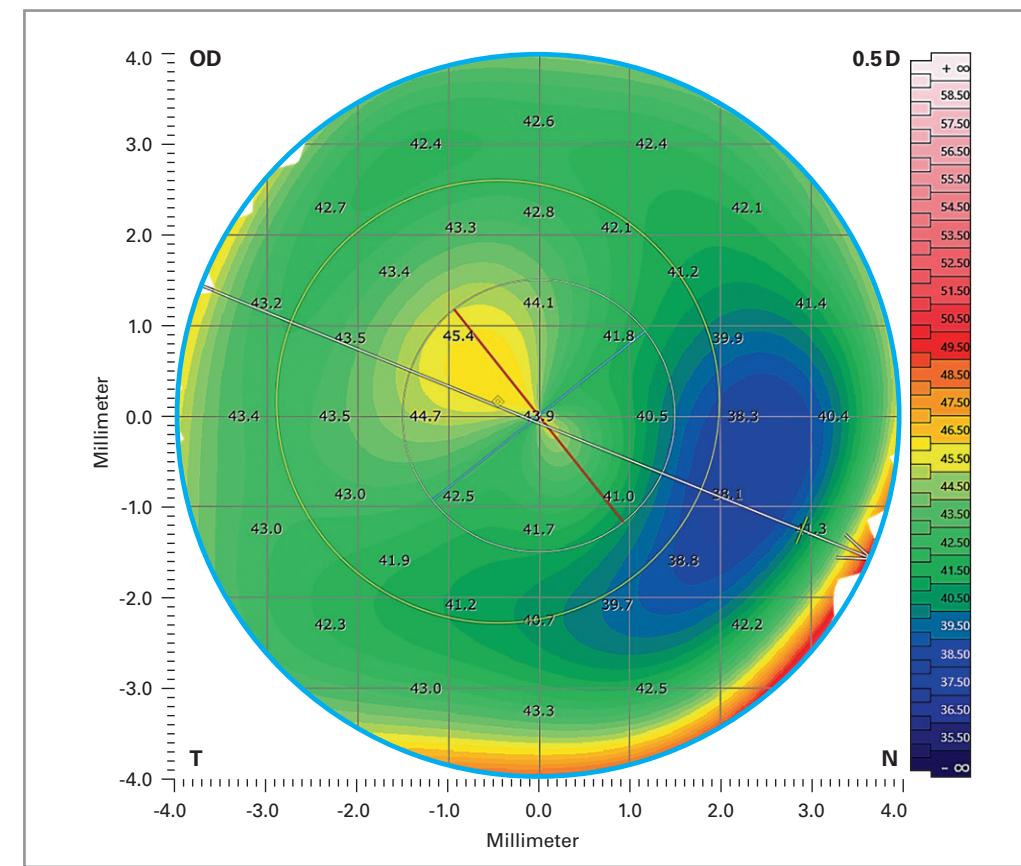

IR camera image

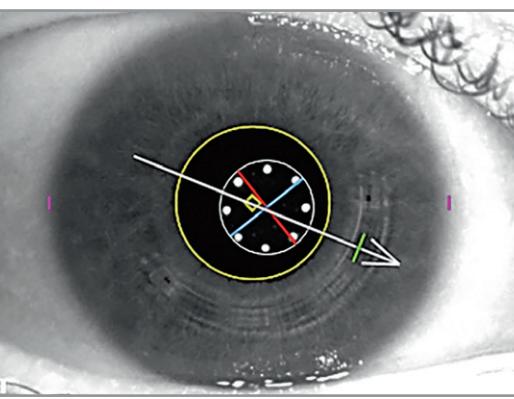

OCT cross-sectional scan

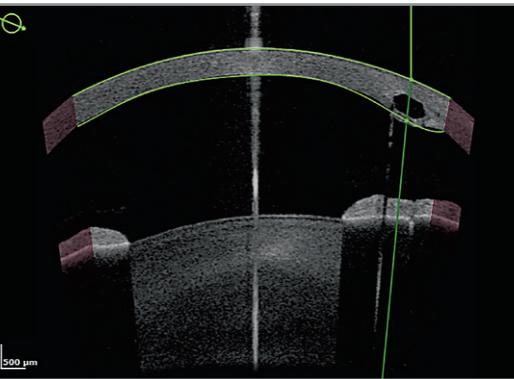


Intrastromal corneal ring segment

Image courtesy: Oliver Findl, MD, MBA, FEBO, Vienna, Austria




Visualize, measure and document progression


Examinations –
Total corneal power map:

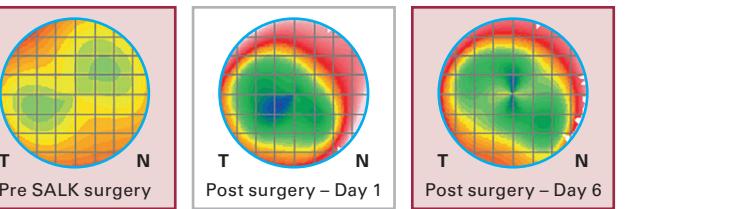

Baseline Pre SALK surgery

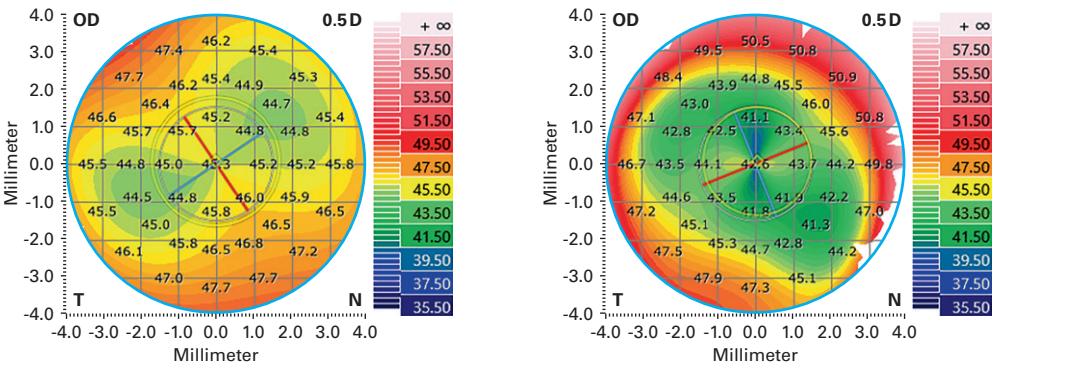


Trend analysis Total corneal power – K (average)

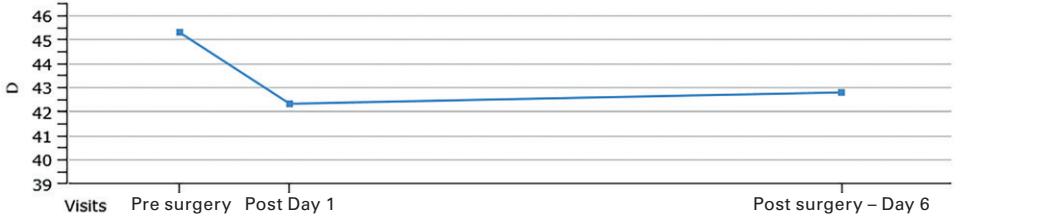


Superficial Anterior Lamellar Keratoplasty – SALK

Image courtesy: Kieren Darcy, MD, BM, MRCS(Eng), CertLRSF, RCOpht, MBA, Bristol, UK

Baseline – Pre SALK surgery

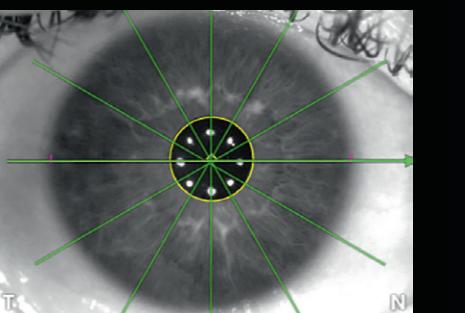

Follow-up – Day 6 –
Post SALK surgery

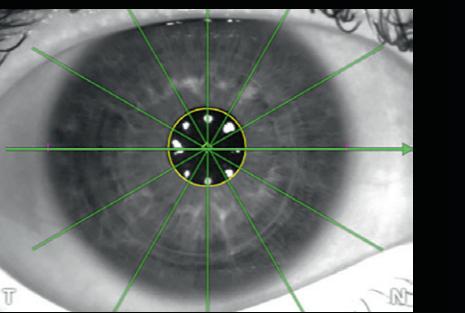


T → N

T → N



Confirm the reliability of your data using the Imaging App if clinically indicated


Image courtesy: Kieren Darcy, MD, BM, MRCS(Eng), CertLRSF, RCOpht, MBA, Bristol, UK

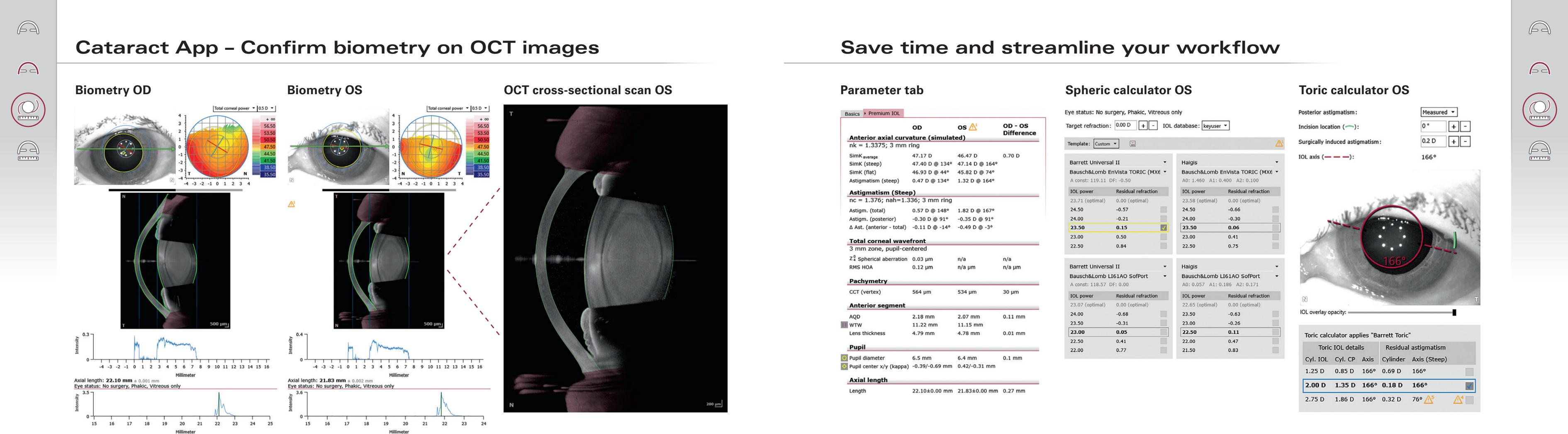

Conveniently compare the biometry data of both eyes. Optimize your clinical workflow using OCT images to visually confirm your biometry data. See what you measure and measure what you see – for more confidence.


Save time and streamline your workflow

Using the dedicated parameter tab, you can easily identify differences between the right and left eye. The integrated spheric and toric IOL calculator adds further convenience to your pre-operative routine and streamlines your workflow.

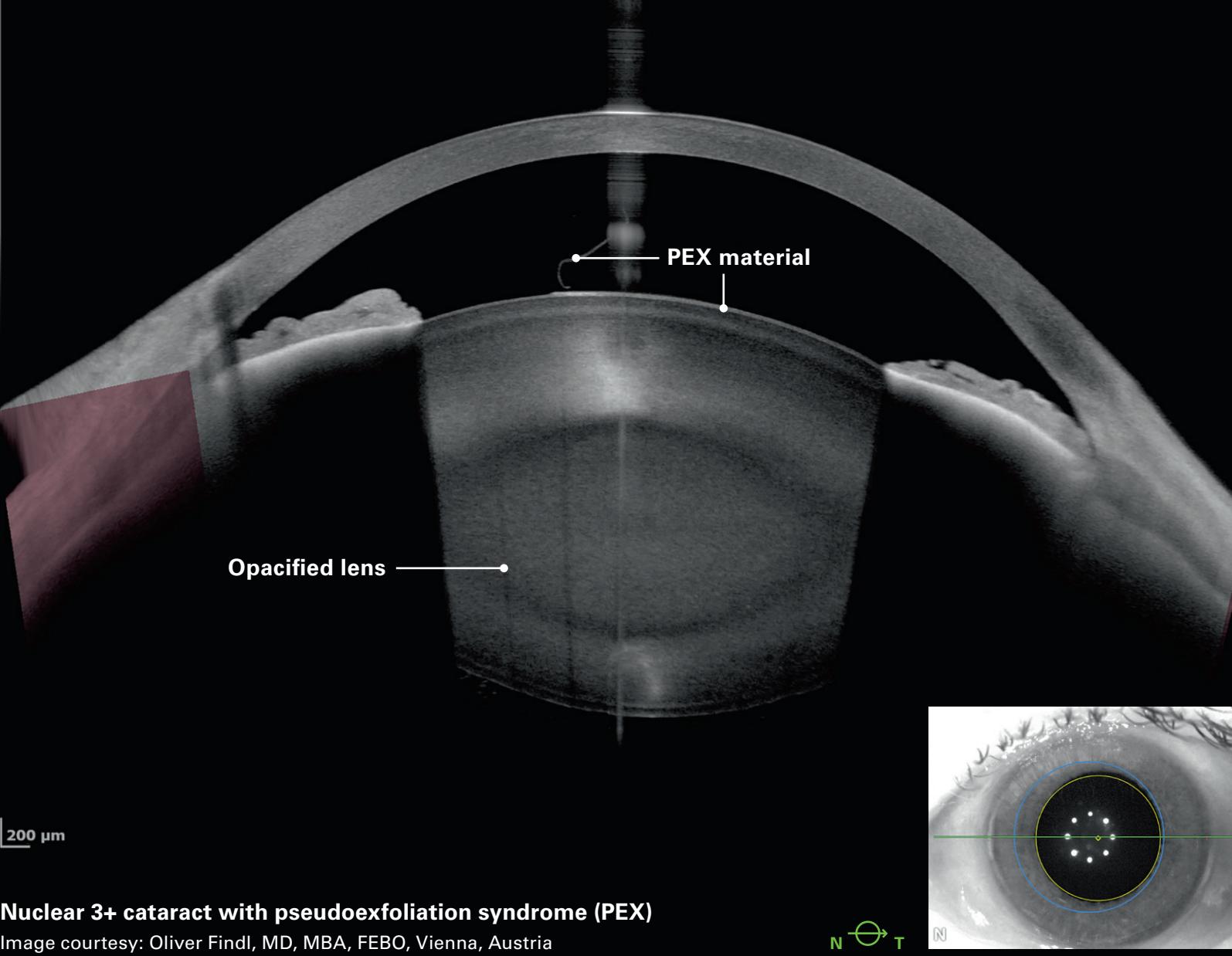

Nuclear 3+ cataract with pseudoexfoliation syndrome (PEX)

Image courtesy: Oliver Findl, MD, MBA, FEBO, Vienna, Austria

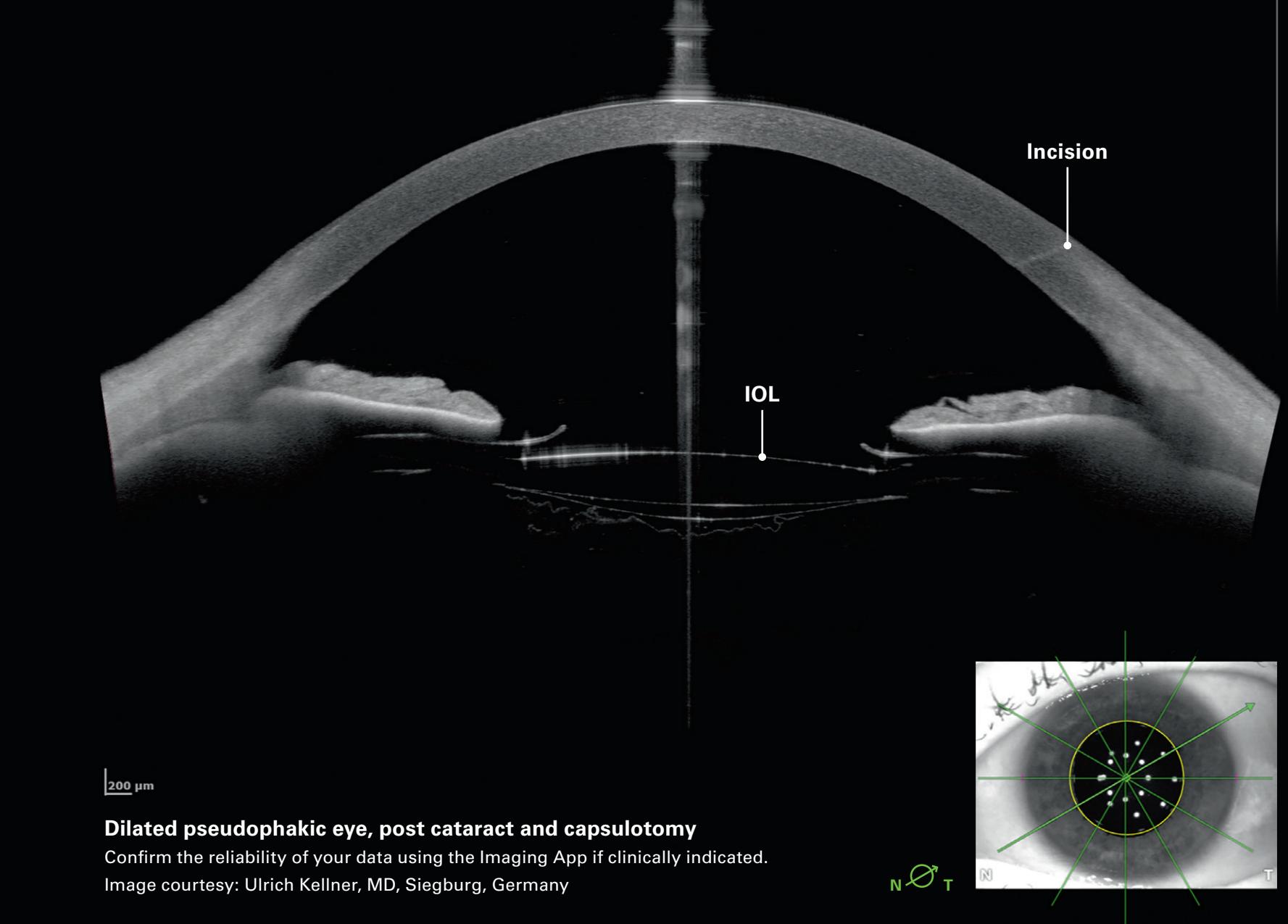


Dilated pseudophakic eye, post cataract and capsulotomy

Confirm the reliability of your data using the Imaging App if clinically indicated.

Image courtesy: Ulrich Kellner, MD, Siegburg, Germany




Metrics App – Measure what you see

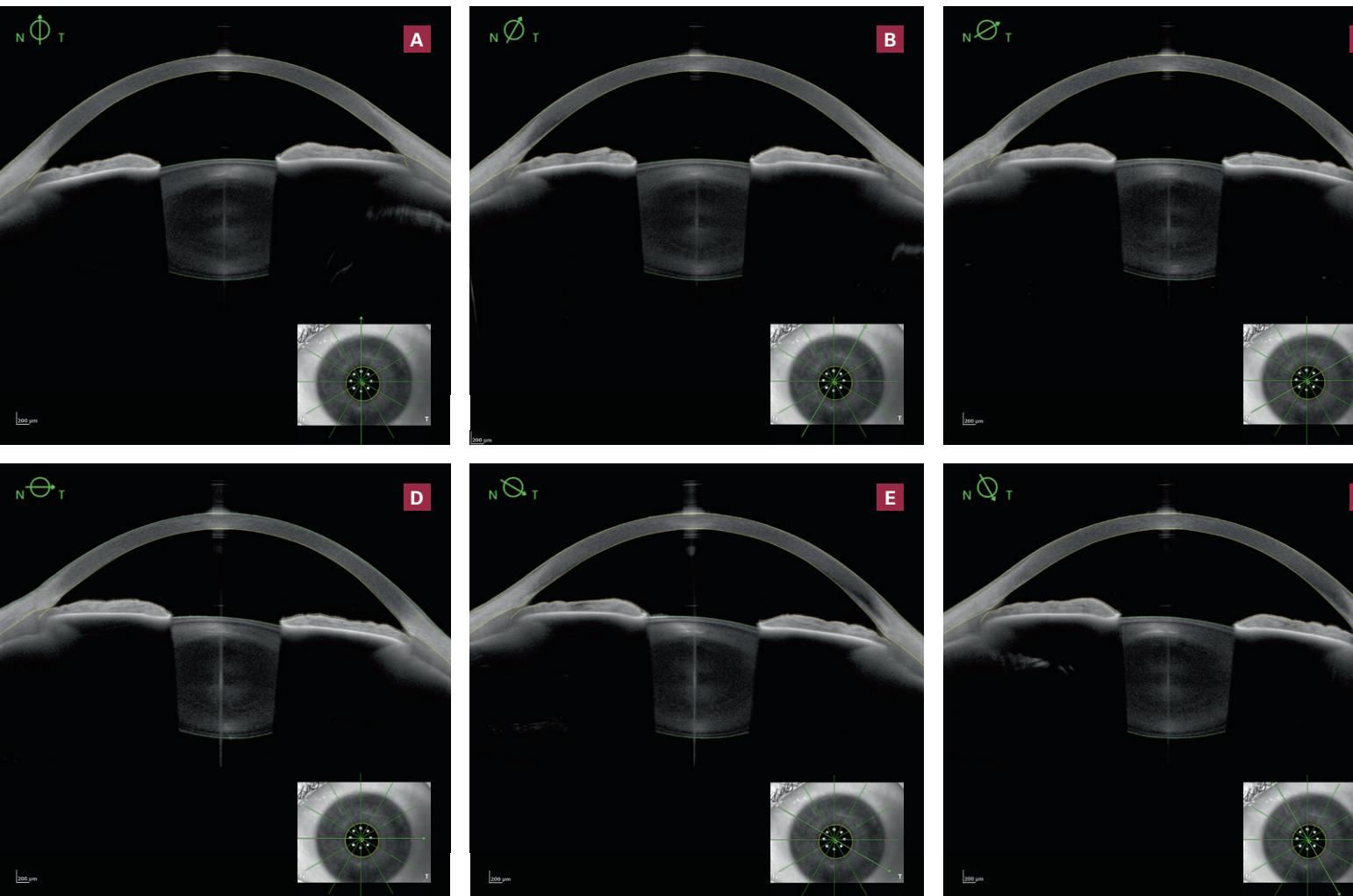


Anterior segment parameters at a glance


ACA 500

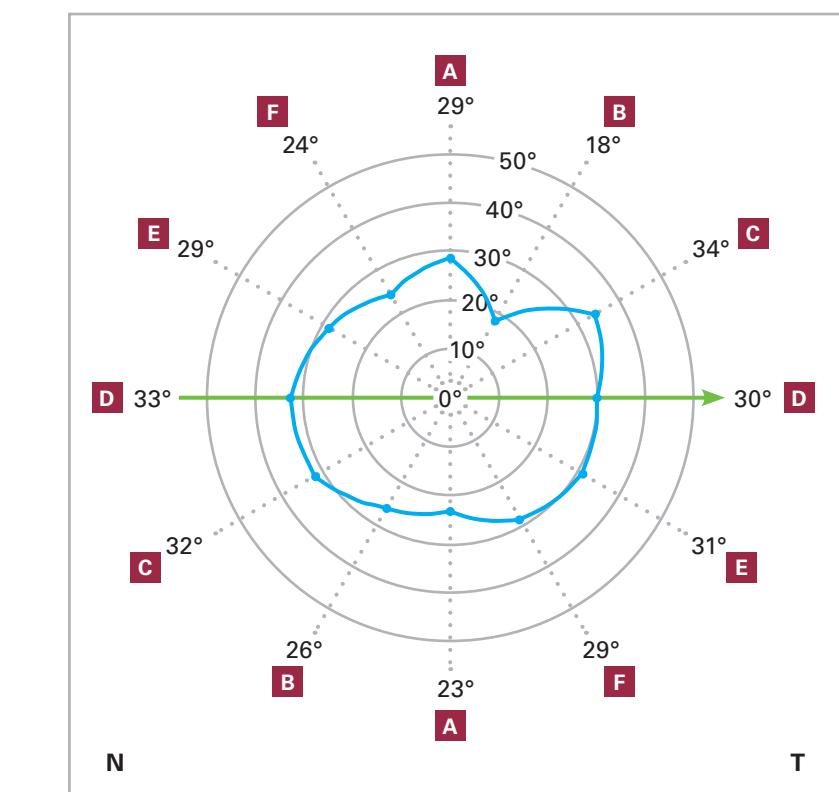

TISA 500

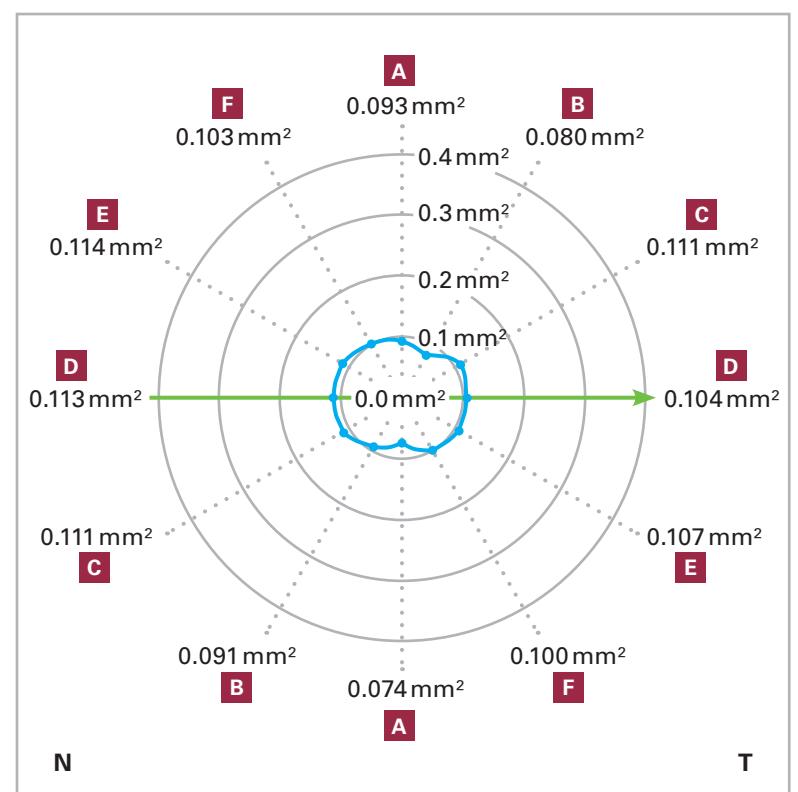

Benefit from anterior segment parameters at a glance using the 360° graphs to optimize your daily clinical routine. The Metrics App features the essential anterior chamber metrics such as aqueous depth, anterior chamber angles (ACA), angle opening distance (AOD), scleral spur angle (SSA), trabecular iris space area (TISA), ACA distance, spur-to-spur distance, central corneal thickness, and white-to-white.

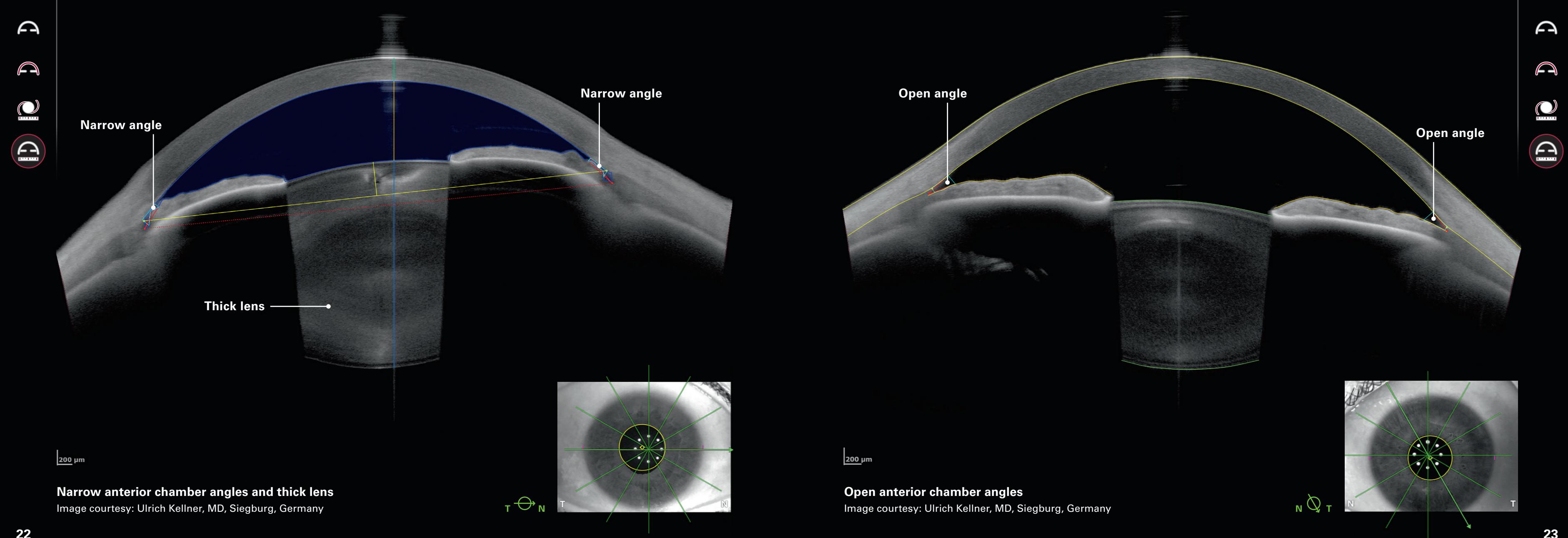



Headquarters

Heidelberg Engineering GmbH · Max-Jarecki-Str. 8 · 69115 Heidelberg · Germany
Tel. +49 6221 64630 · Fax +49 6221 646362

AUS

Heidelberg Engineering Pty Ltd · Suite E5, 63 - 85 Turner St · Port Melbourne VIC 3207
Tel. +61 396 392 125 · Fax +61 396 392 127

CH

Heidelberg Engineering Schweiz GmbH · Schulstrasse 161 · 8105 Regensdorf
Tel. +41 44 8887 020 · Fax +41 44 8887 024

FIN

Heidelberg Engineering GmbH · Luomannotko 6 · 02200 Espoo
Tel. +358 505 226 963

UK

Heidelberg Engineering Ltd. · 55 Marlowes · Hemel Hempstead · Hertfordshire HP1 1LE
Tel. +44 1442 502 330 · Fax +44 1442 242 386

www.HeidelbergEngineering.com